



- IPL
- SIFHR

- RF



Introducing ML600

This is 5 in 1 Multifunctional beauty device. It combines 5 technologies, E-light/IPL & Speed Hair Removal & Nd:yag laser & 808nm Diode Laser & RF. It has 4 handle pieces: Speed Hair Removal handle and 808nm Diode Laser and Nd:yag laser handle and RF handle.



IPL/SHR utilizes the principle of absorption of light to works on hair follicle. The light will be absorbed by Melanin of hair, after the hair follicle be destroyed. Then it will fall off automatically and be removed permanently. Speed Hair Removal also could promote recombination of collagen on the dermis, make the collagen active, and reach function of skin rejuvenation and wrinkle removal.

The ML500 System Highlights



Nd:yag laser erupt high energy in an instant, after absorbing the laser energy, the pigment in the tattoo broke instantly, then becomes tiny pigment group, finally get out from body with metabolism. Meanwhile, it also can do black doll skin rejuvenation treatment, using activated carbons unique characteristics, penetrate into pores, absorbing grease and dirt. Laser light stimulate the regeneration of dermal collagen and the shrink of collagen fibers, to reach a whitening, pigment lightening, pores shrinking and blackheads removal effect, making skin recover young.



808nm Diode Laser damage the chromophores by selectively heating them while leaving surrounding tissue unharmed.

• Permanent hair reduction • For all skin types and skin textures • Including skin rejuvenation treatment



RF can promote subcutaneous collagen contraction, while cooling the skin surface, the dermis layer is heated and the epidermis maintains normal temperature, this time will produce two kinds of reactions: First, the skin dermis thickens, wrinkles follow Lightening or disappearing; second, the remodeling of the subcutaneous collagen, resulting in new collagen, the skin becomes more compact after a treatment.







- Ensure a safe treatment.
- Perfect cooling system: Water, Air.Semiconductor, Display program.
- Easy operation, quick return.
- No down time, no pain.

- Pigmented Lesions
 Skin Remodeling
 Tattoo Removal
 - Acne Removal
 Carbon Peeling
 Hair Removal
 - Vascular Therapy
 Wrinkle Removal



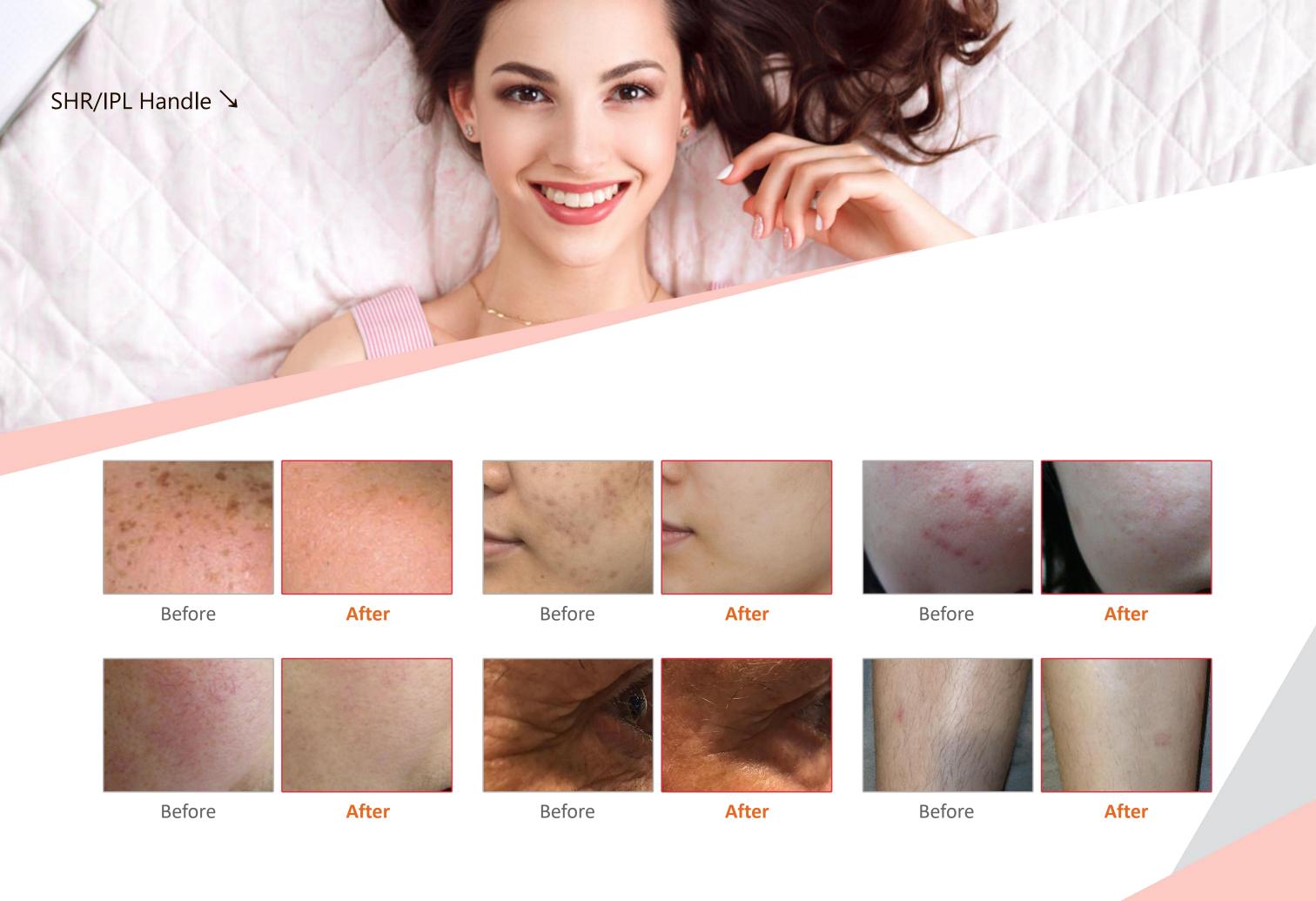

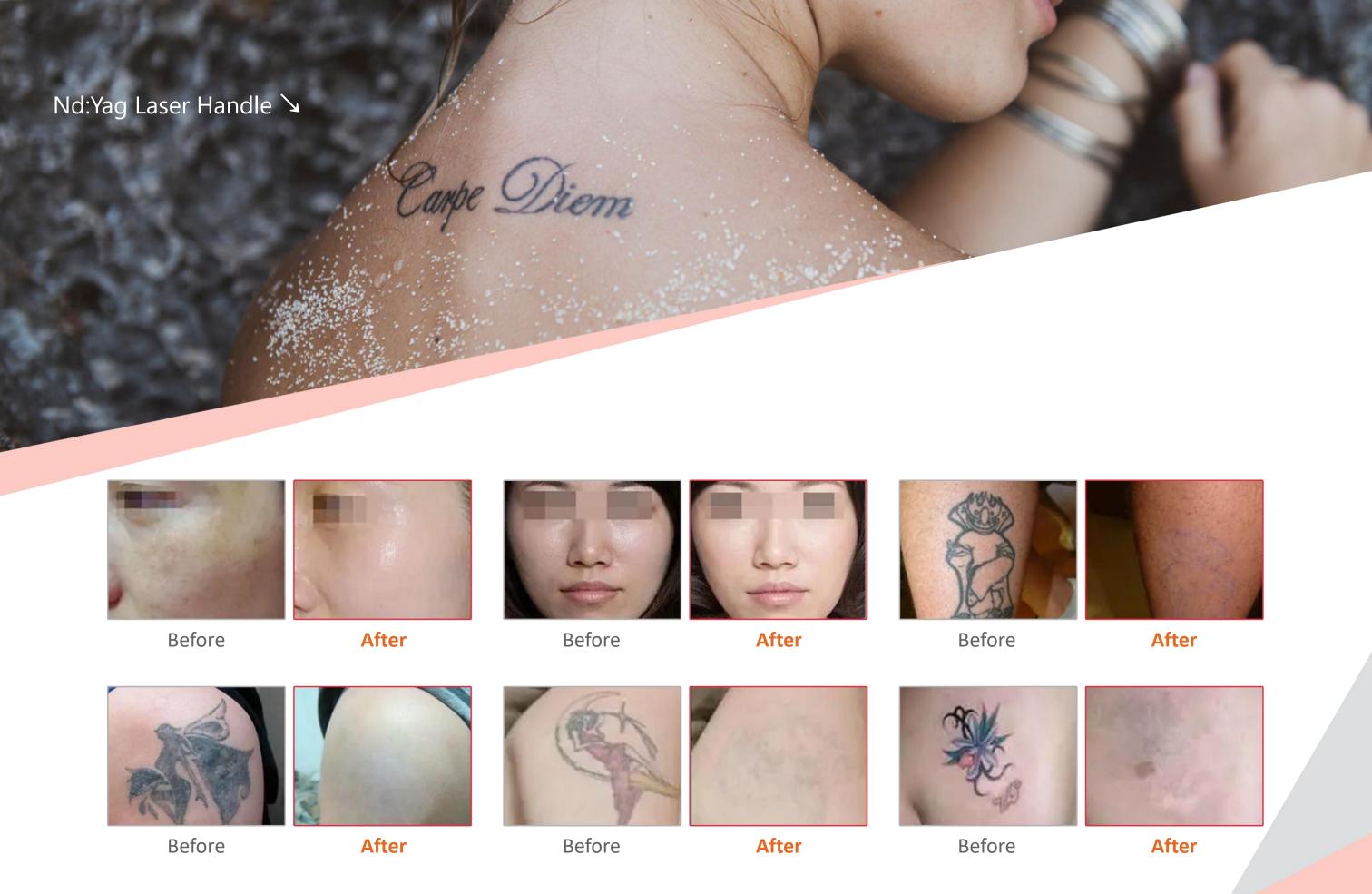















Before

After

Before

After

Before

After













Before

After

Before

After

Before

After















Before

After

Before

After


Before

After







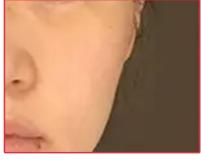





Before

After

Before

After

Before

After

Specification

Model	M600
Language	English, Czech, German21 Languages Optional
Handles	4 handles: SHR, Nd:yag Laser, Diode Laser, RF
SHR/IPL Filters	5 filters(430nm/480nm/530nm/560nm/640nm)
SHR Energy	1-50J/CM2
SHR Frequency	1-10HZ
Nd:yag Laser Energy	510-1000V
Nd:yag Laser Frequency	1-10HZ
Nd:yag Laser Probes	3 probes: 532nm, 1064nm, 1320nm
808nm Laser Device	USA Coherent imported laser bars
808nm Spot Size	14mm*25mm
808nm Energy	1-100W/CM2
808nm Frequency	1-10HZ
RF Probes	20nm, 30nm, 35nm
RF Energy	1-20J
Cooling System	Wind + Water + Semiconductor + HSPC
Output Power	3000W
Machine Size	110*41*59cm
Machine Weight	49KG
Voltage	110V/220V

